

06/2024



MADE IN SWEDEN
www.silentia.com

Product specifications

The Silentia ClearPanels is a transparent, mobile screen designed for easy monitoring and communicating with patients while limiting direct contact. This screen meets strict healthcare hygiene requirements with smooth, hard surfaces that withstand any disinfectants. The wheelbase takes up minimal floor space, while providing maximum stability.

Material

The screens are made of polycarbonate panels and aluminum hinges and frames. The wheelbase consists of a rigid steel base with a plastic cover made of PC/ABS. The wheels are made of polyamide with polyurethane tires. The connections are made of polyamide reinforced with fiberglass. The lightweight construction is user-friendly and optimizes stability and safety.

Dimensions & sizes

Height









1,55 m & 1,85 m

Length

1,50m

Number of panels

3 panels

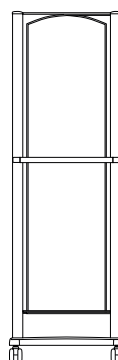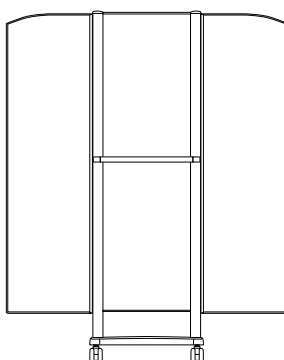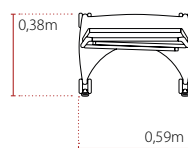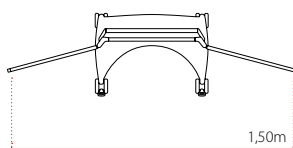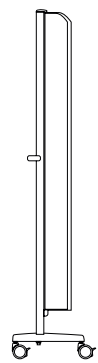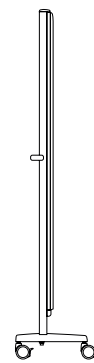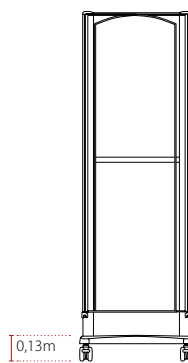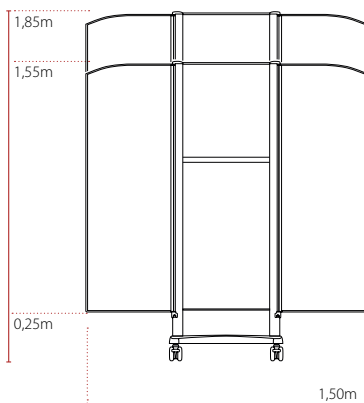

Silentia ClearPanels™

No. of Panels	3	
Length	1,50	
Height	1,55	28 kg
	1,85	31 kg

Product specifications

Care & maintenance

EasyClean™ All products meet strict hygiene requirements and all surfaces are accessible for cleaning. The screens have hard and smooth surfaces that do not trap dust and withstand the majority of detergents that are used in healthcare. The wheels can easily be removed and cleaned in a dishwasher.

Handling/Use

The Clearpanels screen is stable and easy to move, equipped with a single touch-point handle. During use, the wings remain in position with friction brakes. Due to its stable wheelbase, the brakes do not need to be locked at any time.

Sustainability

Silentia focuses on developing long-lasting products that minimize our environmental impact. Since 1990, Silentia has been committed to improving patient care around the world with a safe, hygienic, and resource-efficient privacy system. [LINK TO SUSTAINABILITY PAGE](#)

Recycling & disassembly

The screen is recyclable and all parts can be completely separated and individually delivered to a suitable recycling facility.

Certificates / Listing



EN ISO 9001:2015



EN ISO 14001:2015



Listed as interior solution for the
"WORLD'S GREENEST HOSPITAL"

Tests

Silentia can provide fire tests, durability tests, etc. Please contact the Silentia sales team for specific information.

Warranty

Silentia provides a 2-year warranty. Spare parts are available. Please contact the Silentia sales team for further information about warranty terms and conditions.

Downloads

[LINK TO 2D/3D FILES](#)

[LINK TO SILENTIA MANUAL](#)

[LINK TO CARE & MAINTENANCE VIDEO](#)