


Semifix screen

PRODUCT FACT SHEET

5

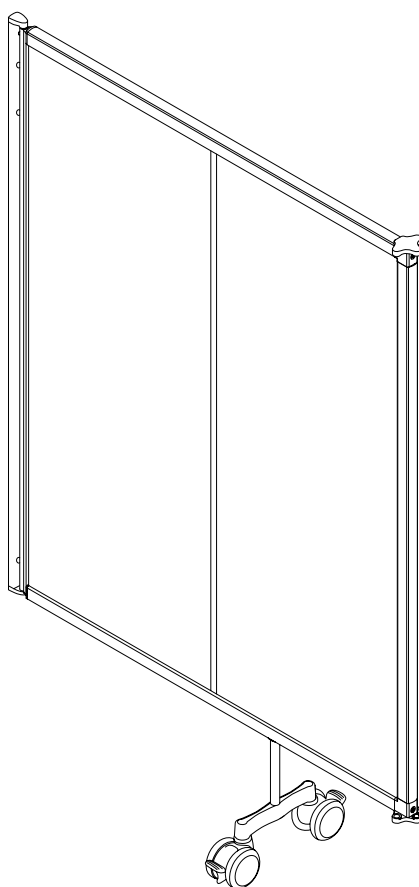

06/2024



MADE IN SWEDEN
www.silentia.com

Product specifications

The Silentia Semifix screen is a large, wall-mounted panel with a low-profile, rotating wheelbase. A three-way connector allows the outside of the Semifix screen to act as a mounting point for up to three additional Semifix or folding screens. With no need for floor attachment, Semifix screens can be used to create private areas within larger open areas. The screens can be attached in succession to form long, sturdy walls or combined with folding screens. The Semifix screen is available in a variety of lengths, heights, and colors.

Material

The screens are made of high-pressure laminate or polycarbonate panels and a aluminum frame. The connections and the dual wheelbase are made of polyamide reinforced with fiberglass. The wheels are made of polyamide with polyurethane tires. The lightweight construction is user-friendly and optimizes stability and safety.









Dimensions & sizes

Height

1,45 m , 1,65 m & 1,85 m

Length

0,75 m to 1,2 m

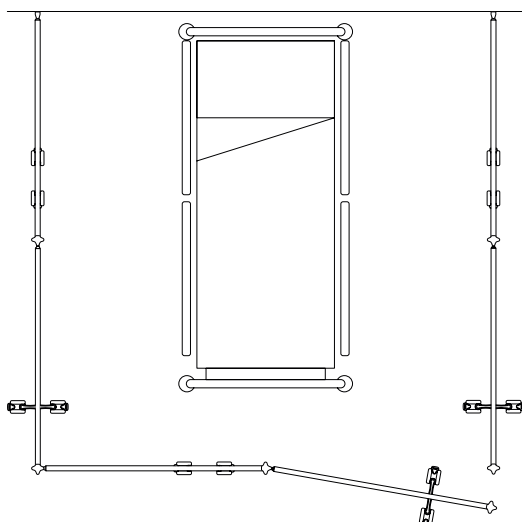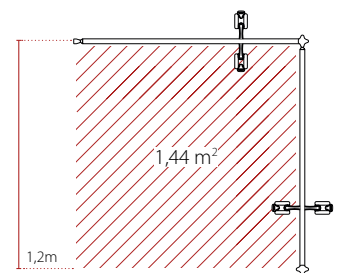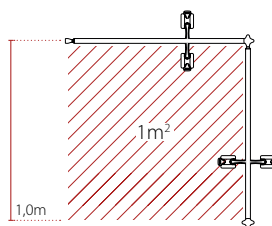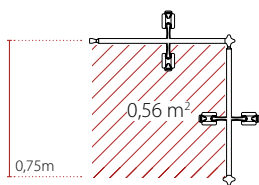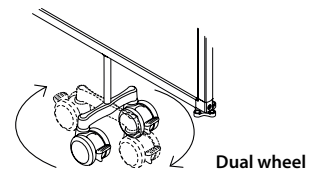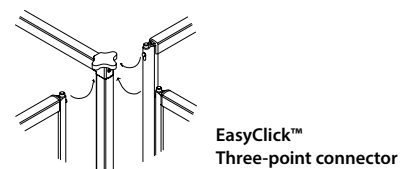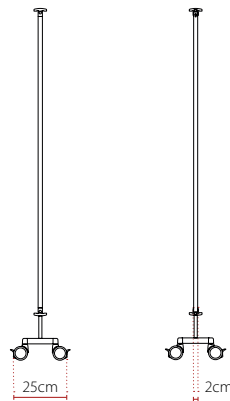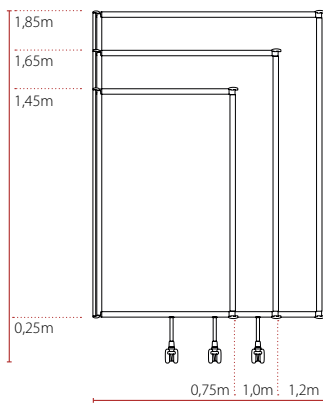

Use of screens example:
6 x Semifix screens (1,2m)
approx. floor area: 5,8m²



Length		0,75	1,00	1,20
Height	1,45	5 kg	6,5 kg	8 kg
	1,65	6,5 kg	8 kg	9,5 kg
	1,85	8,5 kg	10 kg	11,5 kg

Product specifications

Care & maintenance

EasyClean™ All products meet strict hygiene requirements and all surfaces are accessible for cleaning. The screens have hard and smooth surfaces that do not trap dust and withstand the majority of detergents that are used in healthcare. The wheels can easily be removed and cleaned in a dishwasher.

Handling/Use

EasyClick™ Quickly and easily assemble and combine with additional screens. Attach the screen mount to a wall. Connect up to three additional screens with the three-point connector.

The wheelbase of the screen can be rotated 90° for increased stability or extra floor space. The wheels automatically lock in two different positions, perpendicular or parallel with the screen. The wheel suspension compensates for uneven floors.

Sustainability

Silentia focuses on developing long-lasting products that minimize our environmental impact. Since 1990, Silentia has been committed to improving patient care around the world with a safe, hygienic, and resource-efficient privacy system. [LINK TO SUSTAINABILITY PAGE](#)

Recycling & disassembly

The screen is recyclable and all parts can be completely separated and individually delivered to a suitable recycling facility.

Certificates / Listing



EN ISO 9001:2015



EN ISO 14001:2015



Listed as interior solution for the
["WORLD'S GREENEST HOSPITAL"](#)

Tests

Silentia can provide fire tests, durability tests, etc. Please contact the Silentia sales team for specific information.

Warranty

Silentia provides a 2-year warranty. Spare parts are available. Please contact the Silentia sales team for further information about warranty terms and conditions.

Downloads

[LINK TO 2D/3D FILES](#)

[LINK TO SILENTIA MANUAL](#)

[LINK TO CARE & MAINTENANCE VIDEO](#)

Installation

Wall installation



Attaches directly to the wall.


Wall mount directly to wall



Wall mount + spacers




Wall mount on medical rail

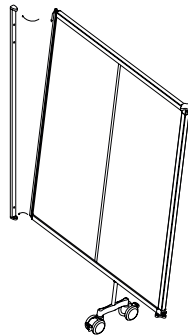




Floor pole installation



For rooms where wall space is limited.
Anchors to the floor.

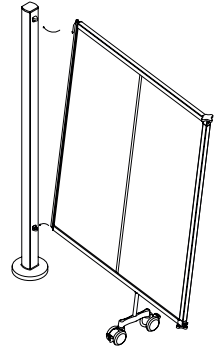

Combination examples

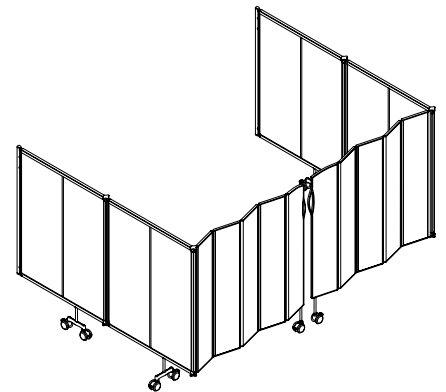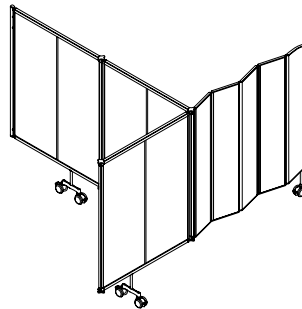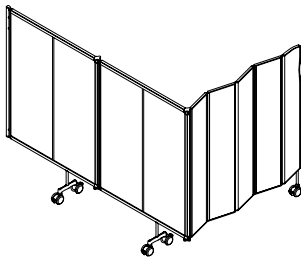

Folding screen



If checked, please fill in PFS 1.

For optimal privacy, combine two Semifix screens with a folding screen.

Folding & fixed screen



If checked, please fill in PFS 1 & 6.


To create two separate private areas, combine three Semifix screens with a folding screen.

Room



If checked, please fill in PFS 1.

To create a room with optimal privacy, combine four Semifix screens with two folding screens.

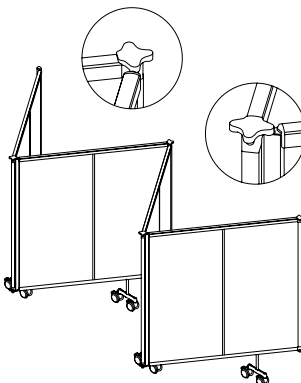

Folding wall




Combination of multiple Semifix screens to build a folding wall.

Design

Color ☐

Laminated panels are available in a variety of colors.

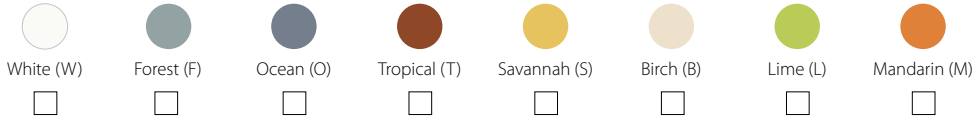

No. of screens: _____


Daylight ☐

Daylight screens feature opaque frosted panels that provide privacy while allowing light into the room.

No. of screens: _____

PhotoPanels™ ☐


PhotoPanels™ create peace of mind and comfort. Suitable for elder care, hospitals, and other healthcare settings. The photos are printed on opaque frosted Daylight panels.

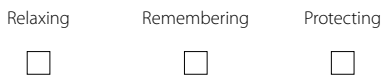

No. of screens: _____ Name of photo: _____

ArtPanels™ ☐

ArtPanels™ is a collection of six custom artist-designed images. The front panel features artwork printed on opaque frosted Daylight panels.

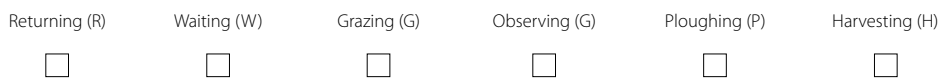

No. of screens: _____

StoryPanels™ ☐

StoryPanels™ provide children with a playful and positive distraction during critical times.

StoryPanels™ Duo combines StoryPanels in pairs with standard laminate panels.



No. of screens: _____

