


Bed end screen - 3 panel

2

PRODUCT FACT SHEET

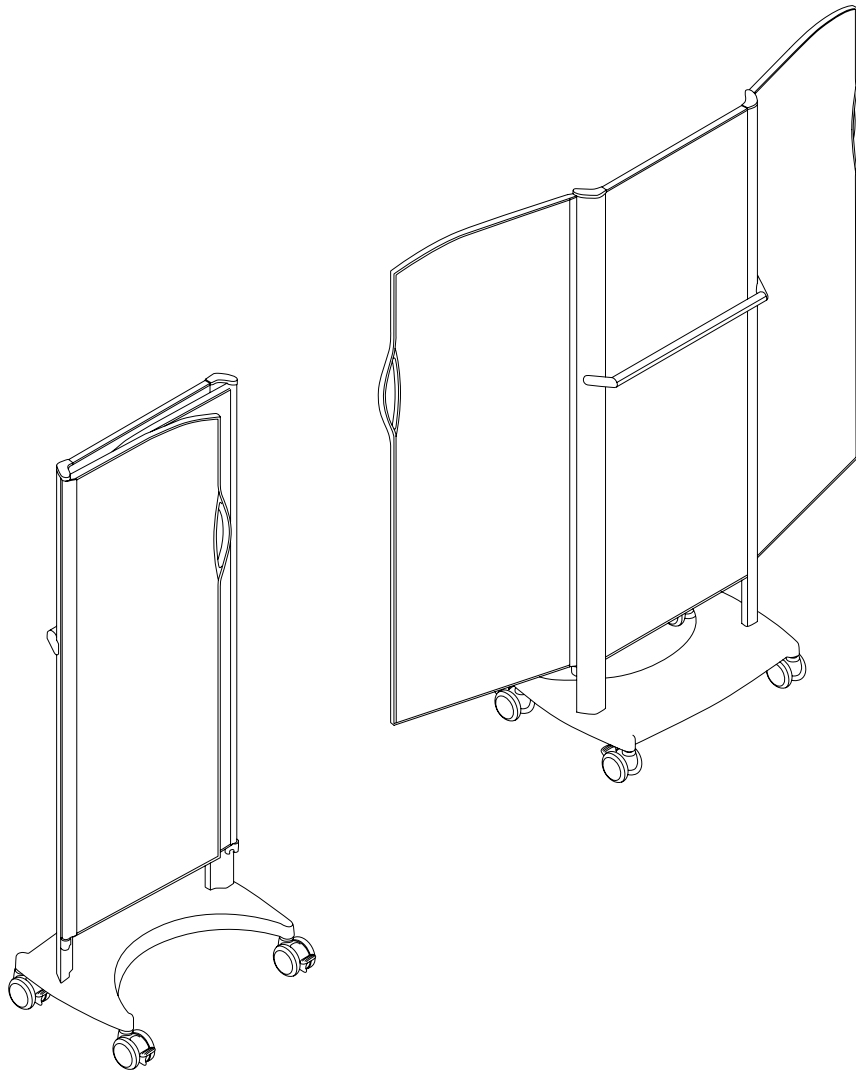

06/2024



MADE IN SWEDEN
www.silentia.com

Product specifications

The Silentia mobile bed end screen is an extremely versatile tri-fold mobile privacy screen. It provides optimal privacy at the end of the bed or whenever a patient needs mobile privacy. The wheelbase takes up minimal floor space, while providing maximum stability. The bed end screen is available in a variety of lengths, heights, and colors.

Material

The screens are made of high-pressure laminate or polycarbonate panels and aluminum hinges and frames. The wheelbase consists of a rigid steel base with a plastic cover made of PC/ABS. The wheels are made of polyamide with polyurethane tires. The connections are made of polyamide reinforced with fiberglass. The lightweight construction is user-friendly and optimizes stability and safety.

Dimensions & sizes

Height









1,55 m & 1,85 m

Length

1,50m

Number of panels

3 panels

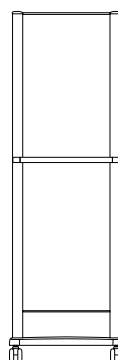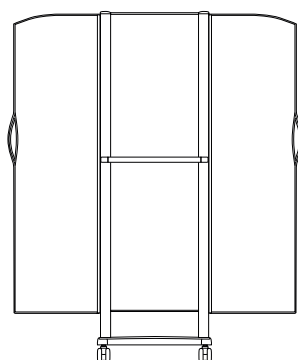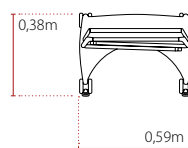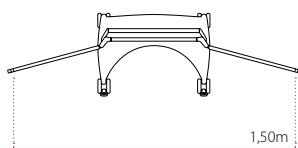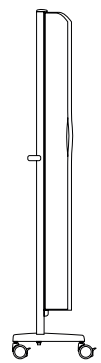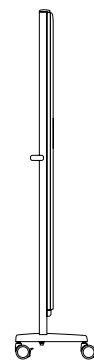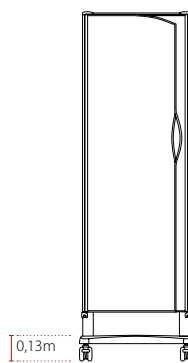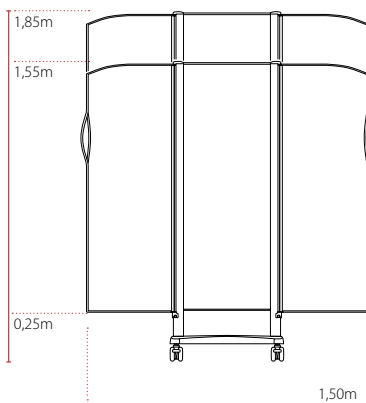

 **Bed end screen**

| | | |
|---------------|------|---------|
| No. of Panels | 3 | |
| Length | 1,50 | |
| Height | 1,55 | 24 kg |
| | 1,85 | 26,8 kg |

Product specifications

Care & maintenance

EasyClean™ All products meet strict hygiene requirements and all surfaces are accessible for cleaning. The screens have hard and smooth surfaces that do not trap dust and withstand the majority of detergents that are used in healthcare. The wheels can easily be removed and cleaned in a dishwasher.

Handling/Use

The bed end screen is stable and easy to move, equipped with a single touch-point handle. The screen can be easily operated from both sides. During use, the wings remain in position with friction brakes. Due to its stable wheelbase, the brakes do not need to be locked at any time.

Sustainability

Silentia focuses on developing long-lasting products that minimize our environmental impact. Since 1990, Silentia has been committed to improving patient care around the world with a safe, hygienic, and resource-efficient privacy system. [LINK TO SUSTAINABILITY PAGE](#)

Recycling & disassembly

The screen is recyclable and all parts can be completely separated and individually delivered to a suitable recycling facility.

Certificates / Listing



EN ISO 9001:2015



EN ISO 14001:2015



Listed as interior solution for the
"WORLD'S GREENEST HOSPITAL"

Tests

Silentia can provide fire tests, durability tests, etc. Please contact the Silentia sales team for specific information.

Warranty

Silentia provides a 2-year warranty. Spare parts are available. Please contact the Silentia sales team for further information about warranty terms and conditions.

Downloads

[LINK TO 2D/3D FILES](#)

[LINK TO SILENTIA MANUAL](#)

[LINK TO CARE & MAINTENANCE VIDEO](#)

Design

Color ☐

Laminated panels are available in a variety of colors.



No. of screens: _____

Daylight ☐

Daylight screens feature opaque frosted panels that provide privacy while allowing light into the room.

No. of screens: _____

PhotoPanels™ ☐


PhotoPanels™ create peace of mind and comfort. Suitable for elder care, hospitals, and other healthcare settings. The photos are printed on opaque frosted Daylight panels.



No. of screens: _____ Name of photo: _____

ArtPanels™ ☐

ArtPanels™ is a collection of six custom artist-designed images. The front panel features artwork printed on opaque frosted Daylight panels.

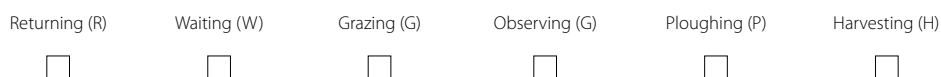


No. of screens: _____

StoryPanels™ ☐

StoryPanels™ provide children with a playful and positive distraction during critical times.

StoryPanels™ Duo combines StoryPanels in pairs with standard laminate panels.

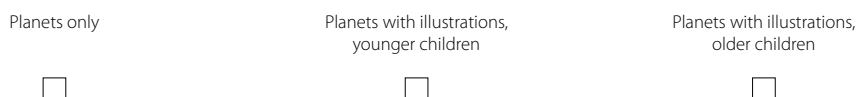

No. of screens: _____

