


# Folding screen

PRODUCT FACT SHEET

1

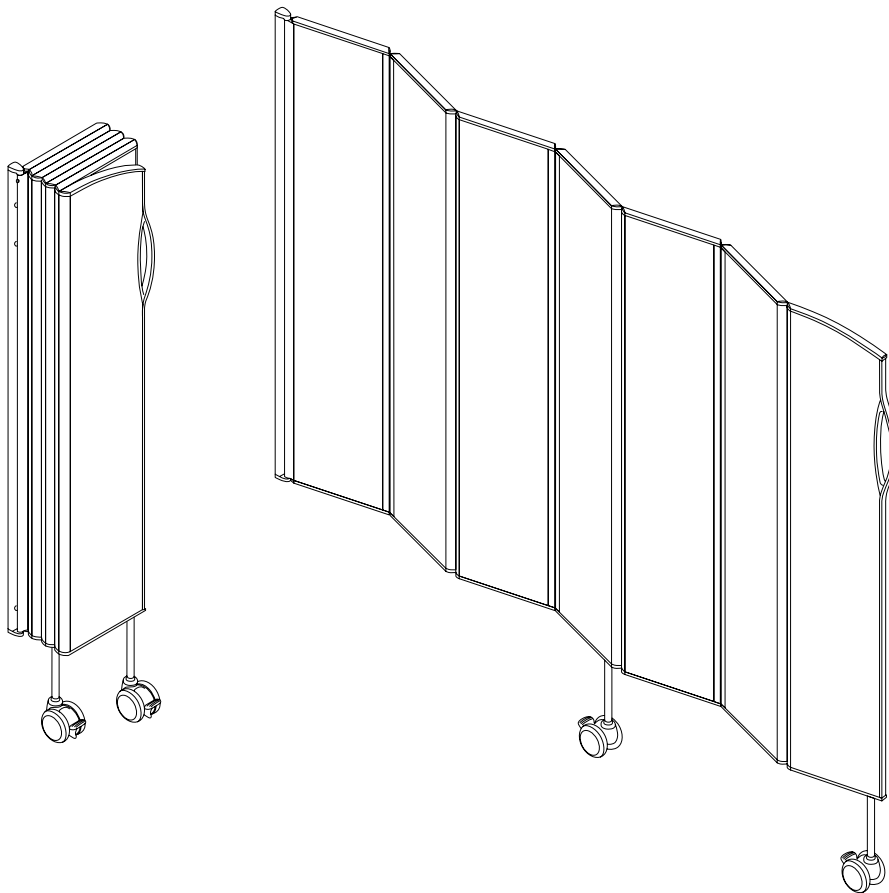

06/2024



MADE IN SWEDEN  
[www.silentia.com](http://www.silentia.com)

## Product specifications

The Silentia folding screen is the essential piece of the Silentia Screen System.

Its accordion design and hard surfaces make it ideal in any care situation, supporting patient privacy and, most importantly, instant cleaning. The folding screen comes in a wide range of lengths, heights, and colors.

### Material

The screens are made of high-pressure laminate or polycarbonate panels and aluminum hinges and frames. The wheels are made of polyamide with polyurethane tires. The lightweight construction is user-friendly and optimizes stability and safety. (5 kg/m<sup>2</sup>)

### Dimensions & sizes

#### Height

1,45m , 1,65m & 1,85m


#### Length

0,75m to 3,75m

#### Number of panels

3 panels to 15 panels

For maximum stability, the longest screens (12-15 panels) are intended for use in an L-shape on a single trolley and U-shape on a double trolley.

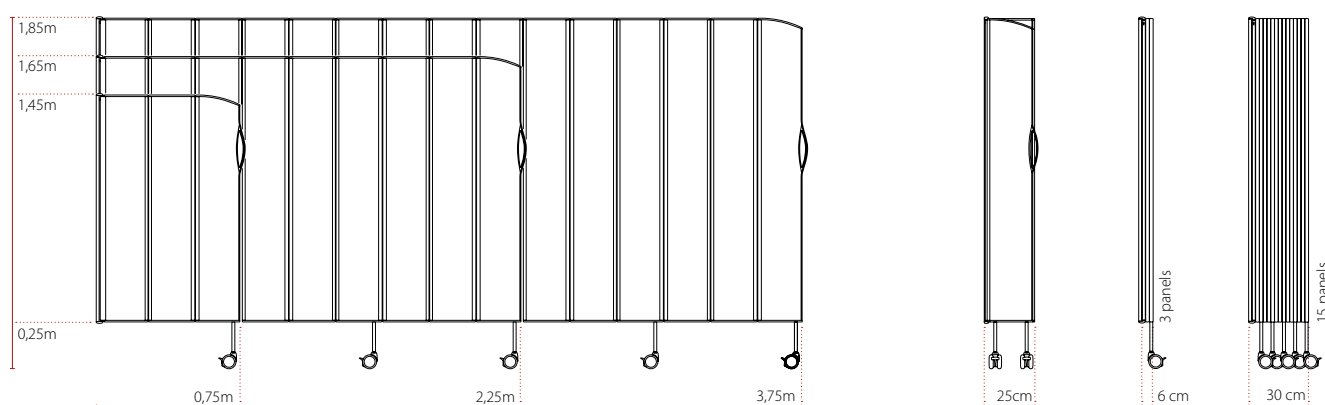

maximum for  
daylight panels  
← →

use as L-shape or  
with dual wheel  
→ ←



No. of Panels										12	13	14	15	
		3	4	5	6	7	8	9	10	11				
Length		0,75	1,00	1,25	1,50	1,75	2,00	2,25	2,50	2,75	3,00	3,25	3,50	3,75
Height	1,45	5 kg	6.5 kg	8 kg	10.5 kg	12 kg	13.5 kg	15.5 kg	17 kg	18.5 kg	20 kg	22 kg	24 kg	25.5 kg
	1,65	6.5 kg	8 kg	9.5 kg	12 kg	14 kg	16 kg	18 kg	19.5 kg	21.5 kg	24 kg	25.5 kg	27.5 kg	29.5 kg
	1,85	8.5 kg	10 kg	11.5 kg	13.5 kg	16 kg	18 kg	20.5 kg	22.5 kg	24.5 kg	27 kg	29 kg	31.5 kg	33.5 kg

# Product specifications

## Care & maintenance

**EasyClean™** All products meet strict hygiene requirements and all surfaces are accessible for cleaning. The screens have hard and smooth surfaces that do not trap dust and withstand the majority of detergents that are used in healthcare. The wheels can easily be removed and cleaned in a dishwasher.

## Handling/Use

**EasyReturn™** The panels of the folding screen are separated by hinges with a patented design that only folds in the intended direction, so that the screen can be folded with just one hand.

**EasyClick™** Regardless of the installation method, a screen can be removed and re-installed on a trolley for mobile use or to any other mount in our system.

## Sustainability

Silentia focuses on developing long-lasting products that minimize our environmental impact. Since 1990, Silentia has been committed to improving patient care around the world with a safe, hygienic, and resource-efficient privacy system. [LINK TO SUSTAINABILITY PAGE](#)

## Recycling & disassembly

The screen is recyclable and all parts can be completely separated and individually delivered to a suitable recycling facility.

## Certificates / Listing



EN ISO 9001:2015



EN ISO 14001:2015



Listed as interior solution for the  
"WORLD'S GREENEST HOSPITAL"

## Tests

Silentia can provide fire tests, durability tests, etc. Please contact the Silentia sales team for specific information.

## Warranty

Silentia provides a 2-year warranty. Spare parts are available. Please contact the Silentia sales team for further information about warranty terms and conditions.

## Downloads

[LINK TO 2D/3D FILES](#)

[LINK TO SILENTIA MANUAL](#)

[LINK TO CARE & MAINTENANCE VIDEO](#)

## Installation

### Wall installation

Choose mounting option below.

Attaches directly to the wall.

Wall mount directly to wall




Wall mount + spacers



Wall mount on medical rail




on

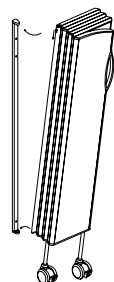


### Floor pole installation

Choose mounting option below.

For rooms where wall space is limited.  
Anchors to the floor.

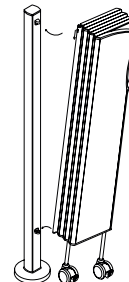


### Combination examples

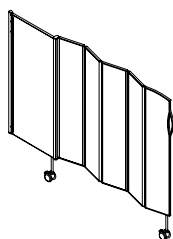


#### Extension screen

If checked, please fill in PFS 4.

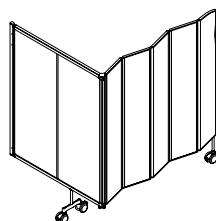

#### Semifix screen

If checked, please fill in PFS 5.

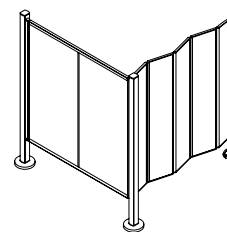


#### Fixed screen

If checked, please fill in PFS 6.




### Mounting point: Right-folding/Left-folding

The mounting point of the folding screen is normally on the left side, which means the screen folds to the right of the wall mount. The screens can also be ordered with the mounting on the right side so the screen folds to the left to avoid interference with a door frame or equipment.

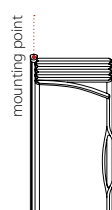

right fold (standard)

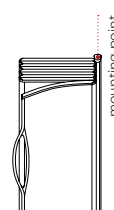


left fold



OR

### Mobile single trolley

Intended for folding screens with 3 to 15 panels.

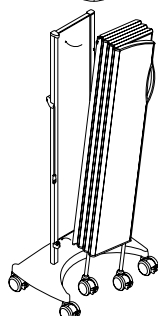


Trolley dimensions:

L: 59 cm / W: 38cm / H: 12,5cm

### Mobile double trolley

A free-standing screen with two folding screens mounted back-to-back.

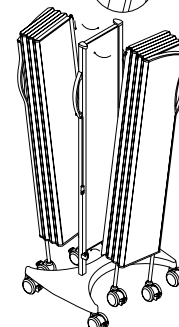

Trolley dimensions:

L: 59 cm / W: 47cm / H: 12,5cm


## Design

### Color ☐

Laminated panels are available in a variety of colors.




No. of screens: \_\_\_\_\_

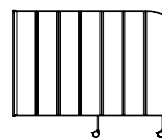

### Daylight ☐

Daylight screens feature opaque frosted panels that provide privacy while allowing light into the room.


No. of screens: \_\_\_\_\_

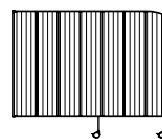

### Daylight Duo ☐




No. of screens: \_\_\_\_\_

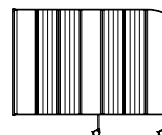

### PhotoPanels™ ☐

PhotoPanels™ create peace of mind and comfort. Suitable for elder care, hospitals, and other healthcare settings. The photos are printed on opaque frosted Daylight panels.

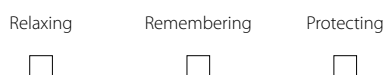


No. of screens: \_\_\_\_\_ Name of photo: \_\_\_\_\_



### ArtPanels™ ☐

ArtPanels™ is a collection of six custom artist-designed images. The front panel features artwork printed on an opaque frosted Daylight panel. The remaining panels are white laminate or Daylight.

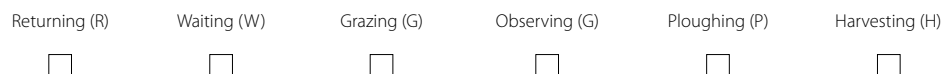


No. of screens: \_\_\_\_\_



### StoryPanels™ ☐

StoryPanels™ provide children with a playful and positive distraction during critical times.

StoryPanels™ Duo combines StoryPanels in pairs with standard laminate panels.

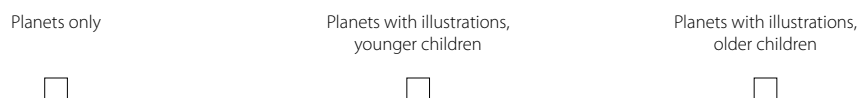

No. of screens: \_\_\_\_\_

